

Automatic Western Blot Processor

WBplex-SA36

Healgen Scientific Limited Liability Company
Address: 3818 Fuqua Street, Houston, TX 77047, USA
Tel: +1713-733-8088
Fax: +1713-733-8848
Website: www.healgen.com
HG-AWB-23001

Automatic Western Blot Processor


WBplex-SA36

Automatic Western Blot Processor is a self-developed instrument. It can automatically complete a series of work such as sealing, incubation, cleaning and recovery of primary antibody and secondary antibody, so as to completely liberate your hands.



Product Characteristics

- The repeatability of injection $CV \leq 3\%$, and the accuracy of injection was $\pm 6\%$.
- Low temperature fluctuations of $\pm 0.3^{\circ}C$.
- Distribution channel: 4-6 channels, with automatic pipeline perfusion and cleaning function.
- Five kinds of alarm prompts: shell opening, manual operation, waste liquid overflow, abnormal liquid absorption and increased membrane strips.

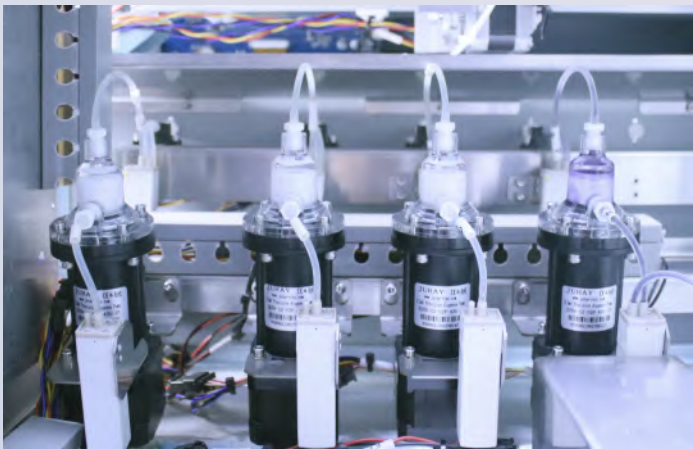

Technical Parameters

Product name	Automatic Western Blot Processor WBplex-SA36
Flux	The instrument has 4 ~ 6 channels and can incubate 36 samples at the same time.
Automation	√
Antibody recovery	√
Temperature controlled system	√
Weight	33kg
Dimension	570mm*510mm*400mm
Precision of liquid's quantities imitted	$CV \leq 3\%$
The edit of program	Flexible sealing, antibody incubation, temperature control and cleaning program settings for faster speed based on actual needs background
Nominal power	Rated power : 120VA

Scope Of Application

- 1、Allergen (Immunoblotting)
- 2、Auto-Antibody (Immunoblotting)
- 3、HIV (Confirmation reagent)
- 4、Research

